Saturn[®]

ST-PS1246

BATHROOM SCALES

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GB

BATHROOM SCALES

Dear Buyer!
We congratulate you on having bought the device under trade name "Saturn". We are sure that our devices will become faithful and reliable assistance in your housekeeping.

Avoid extreme temperature changes. Rapid temperature change (e.g. when the unit is moved from freezing temperature to a warm room) may cause condensation inside the unit and a malfunction when it is switched on. In this case leave the unit at room temperature for at least 1.5 hours before switching it on. If the unit has been in transit, leave it indoors for at least 1.5 hours before starting operation.

This product is specially designed to test body fat/hydration, help keep a close eye on the daily change of your body fat, hydration, and weight, providing information for your reference. Base on the analysis of your height, weight and age etc. it tells your body fat status and gives you an idea of your body condition. Please, read this instruction carefully before use.

■ **User-friendly Tips**
1. During Measuring

A. Always remove your shoes and socks, wipe clean your feet before using.

■ **Performance Index**

Power	Min Weighing Capacity	Max Weighing Capacity	Division	Fat% Range	Fat% Division	Hydration% Range	Hydration% Division	Muscle% Range	Muscle Division
2x3V lithium batteries (CR2032)	2 kg 4.4 lb	150 kg 330 lb	0.1 kg 0.2 lb	4.0% -80.0%	0.1 %	27.5% -86.0%	0.1 %	20.0% -56.0%	0.1 %


■ **Battery**

2x3V lithium batteries (CR2032). Please remove the isolation sheet from the battery contact before use. Observe the polarity.

- B. It is advisable to take measurement at the same time of the day.
- C. Measured result could be misleading after intensive exercise, excessive dieting or under extreme dehydration condition.
- D. Always weigh/use the scale on a hard and flat surface. When you are measuring, please keep still.
- E. The body fat/hydration measured data of the following people may have deviation:
 - Children under 10 or adult over 70 (can only use the device under normal weighing mode)
 - Body building athletes or other occupational athletes.
- F. People with the following conditions should not use the device:
 - Pregnant women
 - Those on dialysis treatment
 - Those who use heart pacemaker or with other implanted medical device.

2. Safety Warning

- A. The Fat% Hydration% Muscle% content evaluated by this device is only for reference (not for medical purpose). If your Fat% Hydration% content is over/below normal level, please consult your doctor for more advice.
- B. For the data accuracy of the scale, please keep it away from strong electro-magnetic field.
- C. It could be slippery if you weigh on a wet surface. Moreover, if the weighing platform surface is glazed, never step on the edge of the scale and do keep balance during your weighing.

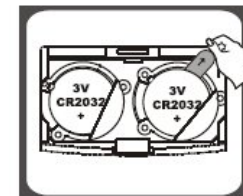
■ To Install Battery


Open the battery case cover at the bottom of scale:

A. Take out the used battery with the aid of a sharp object, as shown in picture.

B. Install new battery by putting one side of the battery beneath the battery baffle and then pressing down the other side.

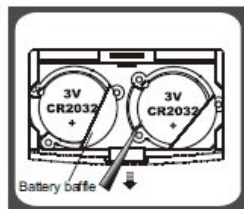■ LCD Display

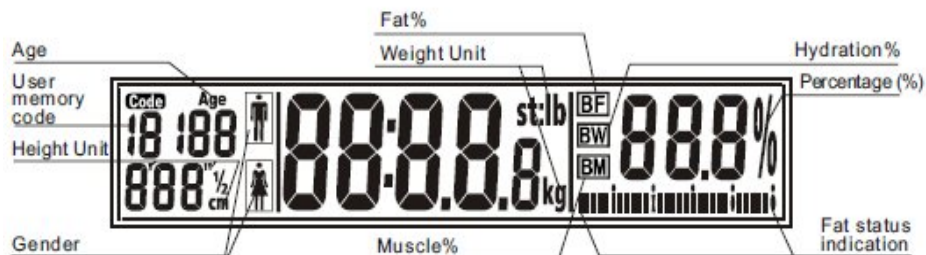
■ To Set Weight Unit

If your scale is available (or applicable) in kg/lb/st, you can set to the one you prefer by pressing the unit conversion button. Press the unit conversion button when scale is power off, LCD shows current weight unit. Press the button again to convert weight unit.

At the bottom of scale you can find a unit conversion button, as shown in picture.

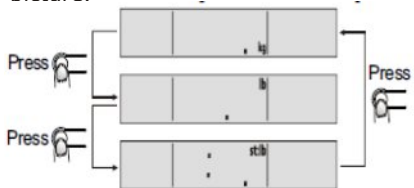After selecting the preferred unit, the scale will be auto-off after 3 seconds. Selected weight unit will be activated when the scale is turned on again.

When setting lb/st as weight unit, the default height unit is feet/inch; When setting kg as weight unit, the default height unit is cm.

■ Key Illustration (Touch button)

SET Turn on the scale and prompt to next parameter. (Please refer to the Setting Personal Parameter overleaf).


▲ Advance the value by one interval each press. Continuous pressing will speed up the advance of value.

▼ Decrease the value by one interval each press. Continuous pressing will speed up the decrease of value.

■ Optional Normal Weighing Mode (bypassing parameter setting)


1. Auto step on weighing

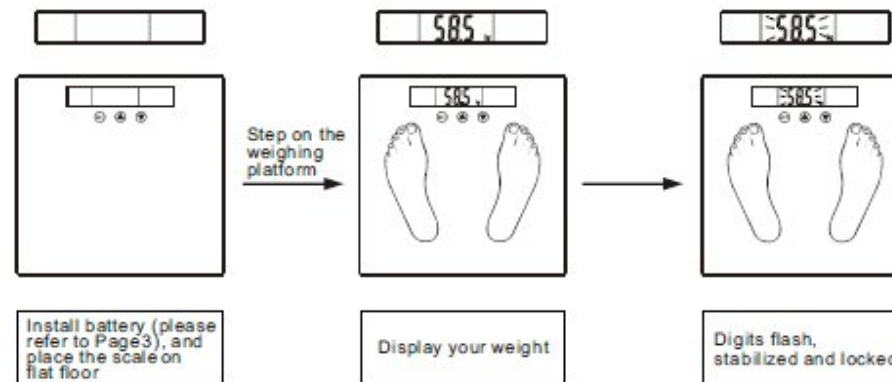2. Consecutive add-on weighing

After digits stabilize, more additional loadings can be conducted. Ensure each add-on loading has to exceed 2kg to get a new reading.

Example: If you want to get a weighing value of your baby, you can weigh yourself first, and then hold your baby together, without stepping off the scale. Then you will get the total weight reading. Subtracting the previous reading of yourself from the new reading, you will get the difference as your baby's weighing value.

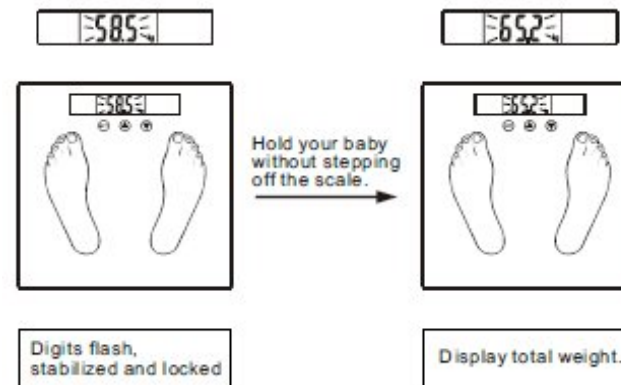■ Fat/Hydration/Muscle Measuring Mode

1. To set personal parameters


① At the first time when you install battery and switch on, all default parameters are set as follow:

Gender	Height	Age
Male	165cm/5'05.0"	25

② Parameters can be set in the following range:

Memory	Gender	Height	Age
1-12	male/female	100-250cm 3'03.5"-8'02.5"	10-100

③ Parameter setting
Example: the 5th user, female, 166cm, 26 years old can set her parameters as follow:

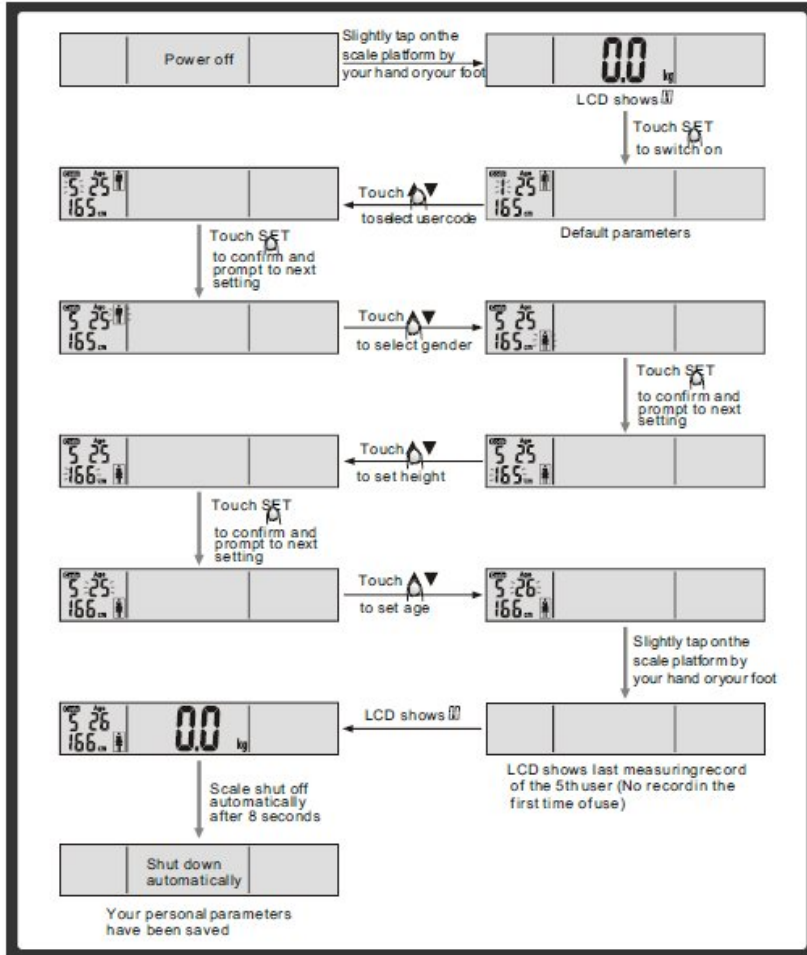
Note:

During parameter setting, you can always tap on the scale platform by your hand or your foot any time to set the scale to zero-point. You can start measuring body fat/hydration based on the parameters just set/updated into the scale. After completing the setting cycle, you can always continue on pressing SET again to go through another setting cycle of all parameters.

2. To Start Measuring


To ensure accuracy, please remove your shoes & socks before stepping on.
Example: the 5th user, female, 166cm, 26 years old can measure as follow:

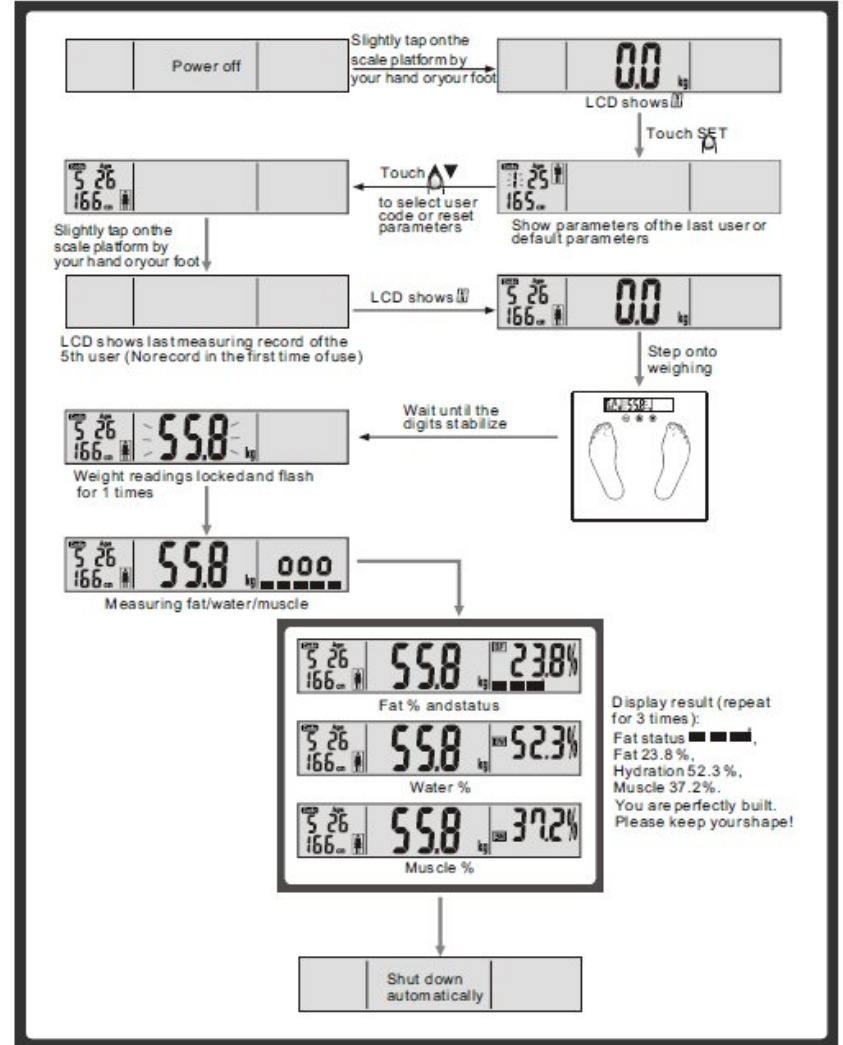■ Warning Indications


Low Battery Indication

The battery power is running low, please replace with a new battery.

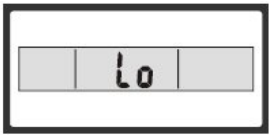
Over-load Indication

The weighing subject on the platform exceeds the maximum capability of scale. Please step off to avoid damage.

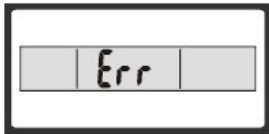
Low Fat% Indication

The Fat% is too low. Please lead a more nutritional diet and take good care of yourself.

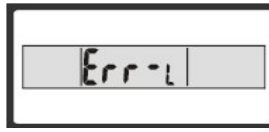
High Fat% Indication

The Fat% is too high. Please watch your diet and do more exercise.

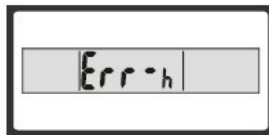Re-test

Error exists, please re-test to get the correct result.

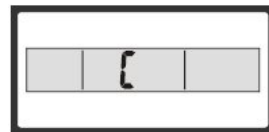■ About Fat

1. About Fat%

Fat% is an index % of fat content in human body.

The fat content in human body is of a certain level, overabundance of fat will bring bad effects on metabolism, thereby, various illnesses may invade.

2. Working Principles

As electrical signal of certain frequency transmits through human body, the impedance of fat is stronger than that of muscle or other human tissues. Through the use of biological resistance analysis, this product will pass a safe frequency of low voltage electrical current through human body, thus, the Fat% will be reflected according to the resistance differences when the electrical signal transmitting through human body.

Cleaning & Maintenance

1. Use a slightly damp soft fabric to wipe clean the scale. DO NOT immerse it into water. Any use of chemical or corrosive detergent is not acceptable.

2. If the plastic components of the scale contacts with grease, spicy sauce, vinegar or other irritating substance, please clean it at once (Please keep the scale away from any fruit juice).

3. Always weight/use the scale on a hard and flat surface. DO NOT use on carpets or mats.

• The actual item may vary from the photographs/illustrations.

3. Fat/Hydration Content Fitness Assessment Chart

After your weighing, you can use the following Assessment Chart for reference:

• The below information is only for reference.

Age	Fat(F)	Hydration(F)	Fat(M)	Hydration(M)	Status	Key Line Indication
≤30	4.0-16.0%	66.0-57.8%	4.0-11.0%	66.0-61.2%	Too lean	
	16.1-20.5%	57.7-54.7%	11.1-15.5%	61.1-58.1%	Lean	
	20.6-25.0%	54.6-51.6%	15.6-20.0%	58.0-55.0%	Normal	
	25.1-30.5%	51.5-47.8%	20.1-24.5%	54.9-51.9%	Fat	
	30.6-60.0%	47.7-27.5%	24.6-60.0%	51.8-27.5%	Too fat	
>30	4.0-20.0%	66.0-55.0%	4.0-15.0%	66.0-58.4%	Too lean	
	20.1-25.0%	54.9-51.6%	15.1-19.5%	58.3-55.3%	Lean	
	25.1-30.0%	51.5-48.1%	19.6-24.0%	55.2-52.3%	Normal	
	30.1-35.0%	48.0-44.7%	24.1-28.5%	52.2-49.2%	Fat	
	35.1-60.0%	44.6-27.5%	28.6-60.0%	49.1-27.5%	Too fat	

Technical Data:

- Ultrafine Design
- Maximum Weight: 150 kg
- Minimum Weight: 2 kg
- Graduation: 100 g
- Fat, Liquid, Muscles Analyzer
- 12 Users Memory
- Age Group: 10-100
- Fat Measurement: 4.0-60.0 %
- Liquid Measurement: 27.5-66.0 %
- Feeding Element Type: 2 Lithium Batteries CR2032

Set

BATHROOM SCALES	1
INSTRUCTION MANUAL WITH WARRANTY BOOK	1
PACKAGE	1

ENVIRONMENT FRIENDLY

DISPOSAL
You can help protect the environment!

Please remember to respect the local regulations: hand in the non-working electrical equipments to an appropriate waste disposal center.

The manufacturer reserves the right to change the specification and design of goods.

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GB

INTERNATIONAL MANUFACTURER'S WARRANTY

The Warranty is provided for the period of 2 years or more if the Law on User Right of the country where the appliance was bought provides for a greater minimum warranty period.

All terms of warranty comply with to the Law on Protection of the Consumer Rights and are governed by the legislation of the country where the product was bought.

The warranty and free of charge repair are provided in any country where the product is delivered to by Saturn Home Appliances or authorized representative thereof, and where no restriction of import or other legal provisions interfere with rendering of warranty services and free of charge repair.

Cases uncovered by the warranty and free of charge repair:

1. Guarantee coupon is filled in improperly.
2. The product has become inoperable because of the Buyer's non-observance of the service regulations indicated in the instruction.
3. The product was used for professional, commercial or industrial purposes (except for the models, specially intended for this purpose that is indicated in the instruction).
4. The product has external mechanical damage or damage caused by penetration of the liquid inside, dust, insects and other foreign objects.
5. The product has damages caused by non-observance of the rules of power supply from batteries, mains or accumulators.
6. The product was exposed to opening, repair or modification of design by the persons, unauthorized to repair; independent cleaning of internal mechanisms etc. was made.
7. The product has natural wear of parts with limited service life, expendables etc.
8. The product has deposits of scale inside or outside of fire-bar elements, irrespective of used water quality.
9. The product has damages caused by effect of high (low) temperatures or fire on non-heat resistant parts of the product.
10. The product has damages of accessories and nozzles which are included into the complete set of product shipment (filters, grids, bags, flasks, cups, covers, knives, sealing rings, graters, disks, plates, tubes, hoses, brushes, and also power cords, headphone cords etc.).

On the warranty service and repair, please contact the authorized service centers of "Saturn Home Appliances".

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GB

WARRANTY COUPON

**When purchasing the product,
please require its checking before you;**

BE SURE

**that the goods sold to you, are functional and complete
and that the warranty coupon is filled in correctly.**

This warranty coupon confirms the absence of any defects in the product you purchased and provides for free of charge repair of the product failed through the fault of the manufacturer throughout the period of warranty service and free of charge repair.

Failing the presentation of this coupon, in case of its improper filling in, infringement of factory seals (if any), and also in cases indicated in the warranty obligations, the claims are not accepted, and no warranty service and free of charge repair is made!

The warranty coupon is valid only in the original copy with the stamp of trading organization, signature of the seller, date of sale and signature of the buyer.

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